


ORTHOPLUS SUPPORTER

QUICKCHART

REAR BUTTON' COLORS



STATIONARY MODE = BLUE
Locks in the height chosen


+ up  - down


SETTING SWING LIMITS

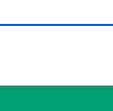
1 LONG CLICK → **1** → **2** → **3** →

Set the **upper** swing limit Set the **lower** swing limit Back in **STATIONARY** mode


↓ ↓

 Position the device where you would like to set the limit with the control pad


 To **SET**, **1 LONG CLICK**
Long vibration = **OK**

 To skip **WITHOUT** setting, **1 SHORT CLICK**

FREE MODE = GREEN
Accompanies the movements



+ **+100g / clic**  - **-100g / clic**

SETTING AN IDEAL COMPENSATION FORCE

 **1 LONG CLICK** on the rear button
Short vibration = **OK**

SLEEP MODE = BLINKING WHITE
2 SHORT CLICKS on the rear button to go back in SLEEP mode

CHANGE THE MODE

 **1 LONG CLICK** on both buttons of the control pad
OR
 **1 SHORT CLIC** on the rear button
Short vibration = **OK**

PARKING POSITION ON WHEELCHAIR

In **STATIONARY** mode, pull down the ORTHOPUS Supporter on the armrest: the device stays blocked in place.



Rear button
blinking blue

Parking position activated = OK

To go out of the Parking position, click on the **+** button of the control pad.

SAFETY IN FREE MODE



Rear button
blinking green

If the device is switched in FREE mode without arm in the arm brace, the ORTHOPUS Supporter puts itself in safety.

WARNING LIGHTS



Unplug the device and wait a few minutes before using it again. If the Rear button stays red, please contact your ORTHOPUS Supporter Country Representative.

RECOMMENDATIONS

i While starting the **ORTHOPUS Supporter**, don't touch the control pad and don't install an arm in the arm brace.

i Prior to removing the user's arm from the support, please **always ensure** that the ORTHOPUS Supporter is in STATIONARY or SLEEP mode.

i It is recommended that users **become acclimated to the device gradually** and be monitored by a qualified healthcare professional.

i Only an **individual trained specifically** is authorized to install the ORTHOPUS Supporter.

To take full advantage of all that the ORTHOPUS Supporter has to offer, **please make sure:**

✓ In horizontal position, the arm isn't leaning forwards or backwards.

✓ The arm is always in contact with the elbow support. If necessary, set an upper swing limit.

✓ On the wheelchair, cables and box interface are protected for the rain.

WARNINGS

⊖ The ORTHOPUS Supporter is **designed exclusively to support the arm**: it is not to be used as an aid when standing up or sitting down, or for any other purpose.

⊖ In the event of faulty or damaged on electronics parts, **please refrain from using the device**. Please contact your ORTHOPUS Supporter Country Representative for any maintenance issues.

⚠ Do not place the device in direct sunlight or directly near a source of heat for an extended period of time.

⚠ The device is water-resistant, but **not waterproof**. Do not expose it to heavy rainfall or significant ambient air humidity.



If the ORTHOPUS Supporter isn't mounted, **it must remain in its packaging and box at all times**. The box is to be kept should the device need to be returned.



Please refer to the ORTHOPUS Supporter User manual on orthopus.com/en/documentation



ORTHOPUS Supporter Quickchart
v1.1 - 2023/03/10